

Minimizing Pressure In-homogeneities for Large Samples in High Pressure Neutron Scattering Measurements

Juscelino B. Leão

NIST Center for Neutron Research

New Directions for High-Pressure Neutron Workshop

June 3–5, 2013

Oak Ridge, TN USA

Minimizing Pressure In-homogeneities for Large Samples in High Pressure Neutron Scattering Measurements


- Reasoning
- Apparatus
- Results
- Technique

Reasoning

- A. Freezing (V_{constant})
- B. Cooling to Freezing Point (P_{constant})
- C. Freezing and Cooling (VP_{constant})


W.F. Sherman, A.A. Stadtmuller. *Exp. Tech. in H. Pressure Research*

- i. Pressurize above P_xT curve
- ii. Cool slowly under P_{constant}
down to the freezing point
- iii. Hope for the best down to base

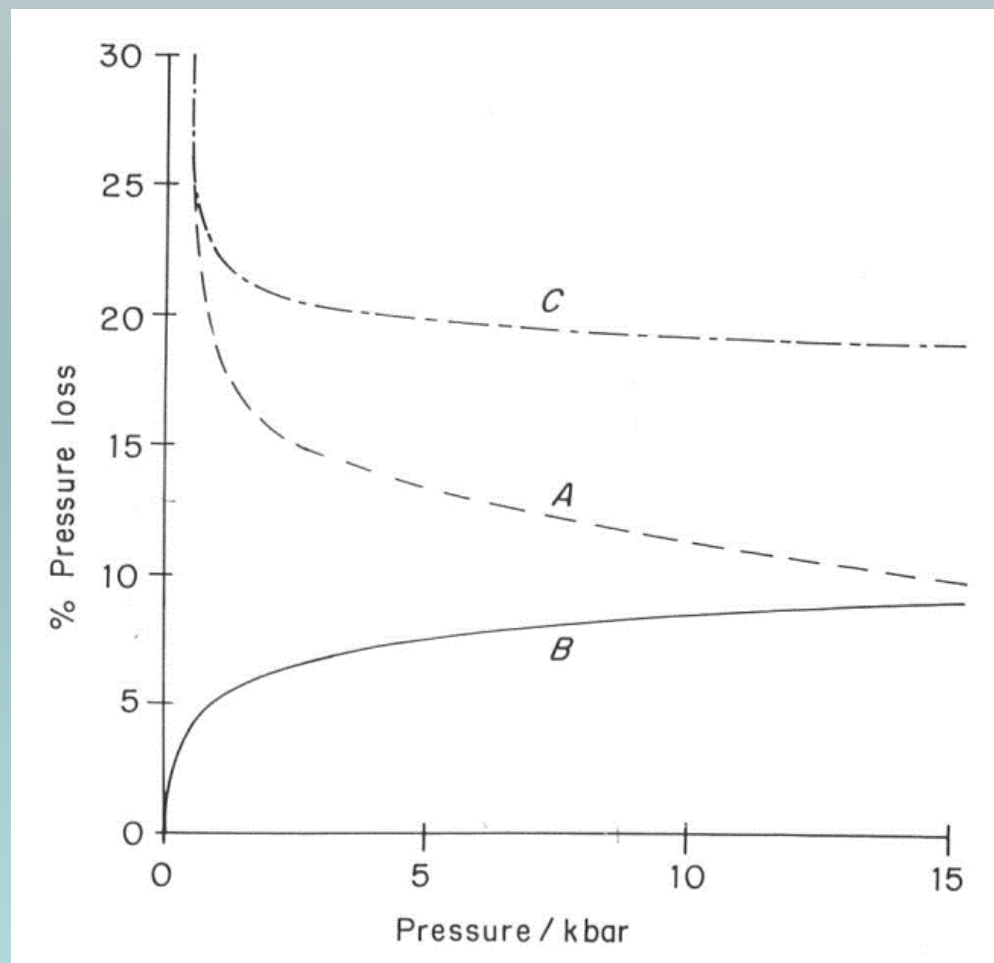

W.F. Sherman, A.A. Stadtmuller. *Exp. Tech. in H. Pressure Research*.

Apparatus

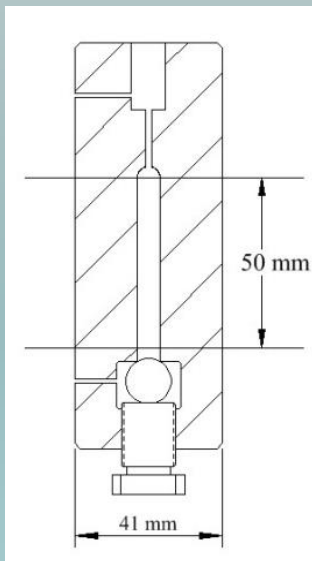

$P_{\text{max}} = 7.0 \text{ kbar}$

Working Pressure = 6.5 kbar

Al 7075-T6 Construction

1.5 cm³ sample volume

69% Neutron transmission at 2Å



Harwood Eng., Inc. 2-Stage Intensifier

Results

- i. Pressurize above P_xT curve
- ii. Cool slowly under V_{P_{constant}}
down to the freezing point
- iii. Hope for the best down to base


$$P = \left(\frac{\beta_o}{\beta'} \right) \left[\left(\frac{r}{r_o} \right)^{-\beta'} - 1 \right]$$

From 1-D Analog to the Murnaghan Equation

$$\beta' = 10.8(9)$$


$$\beta_o^{-1} = - \left(\frac{d \ln r}{dP} \right)_{P=0} = 373^{-1} \text{ kbar}$$

Hanfland, Beister, Syassen. *Phys. Rev. B* **39**, 1989

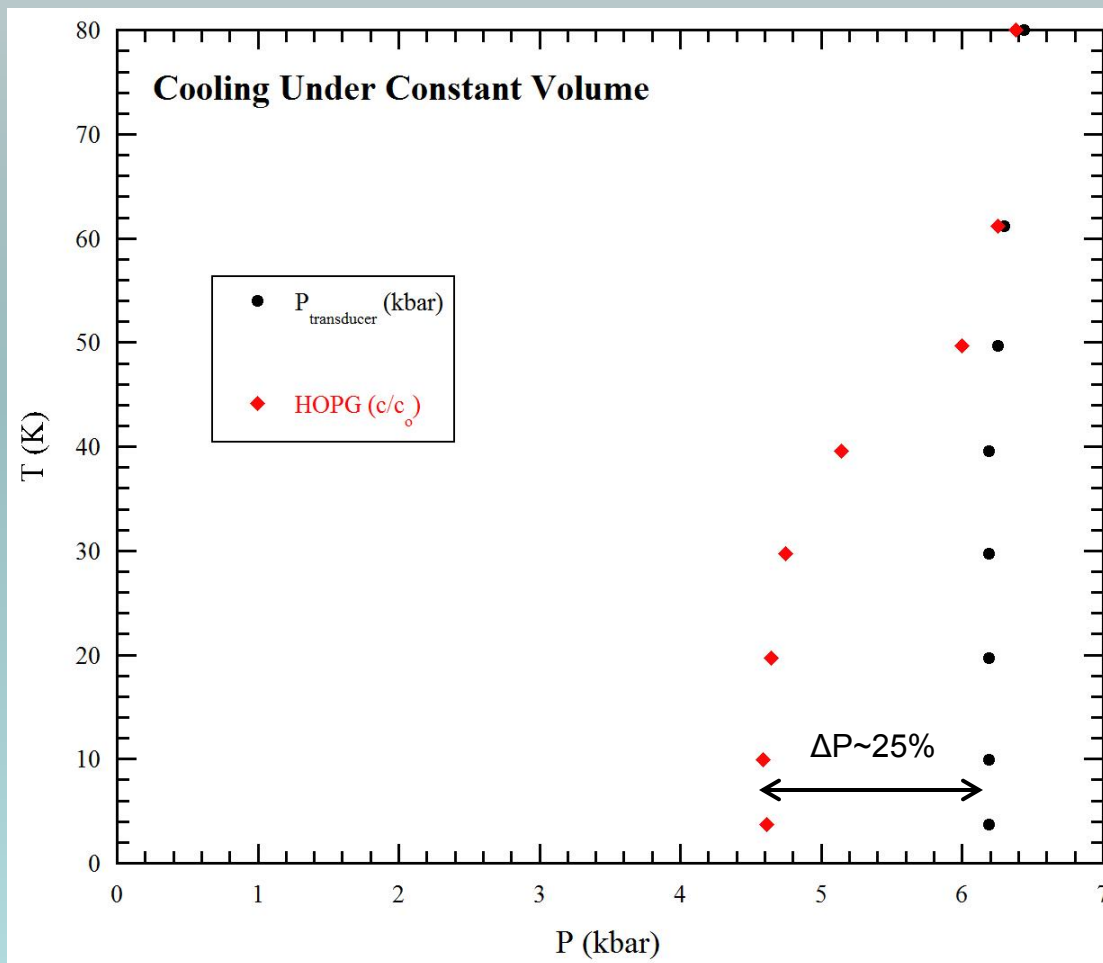

Results


Freezing and Cooling (VP_{constant})
for neutron elastic measurements of
HOPG (c/c_0)

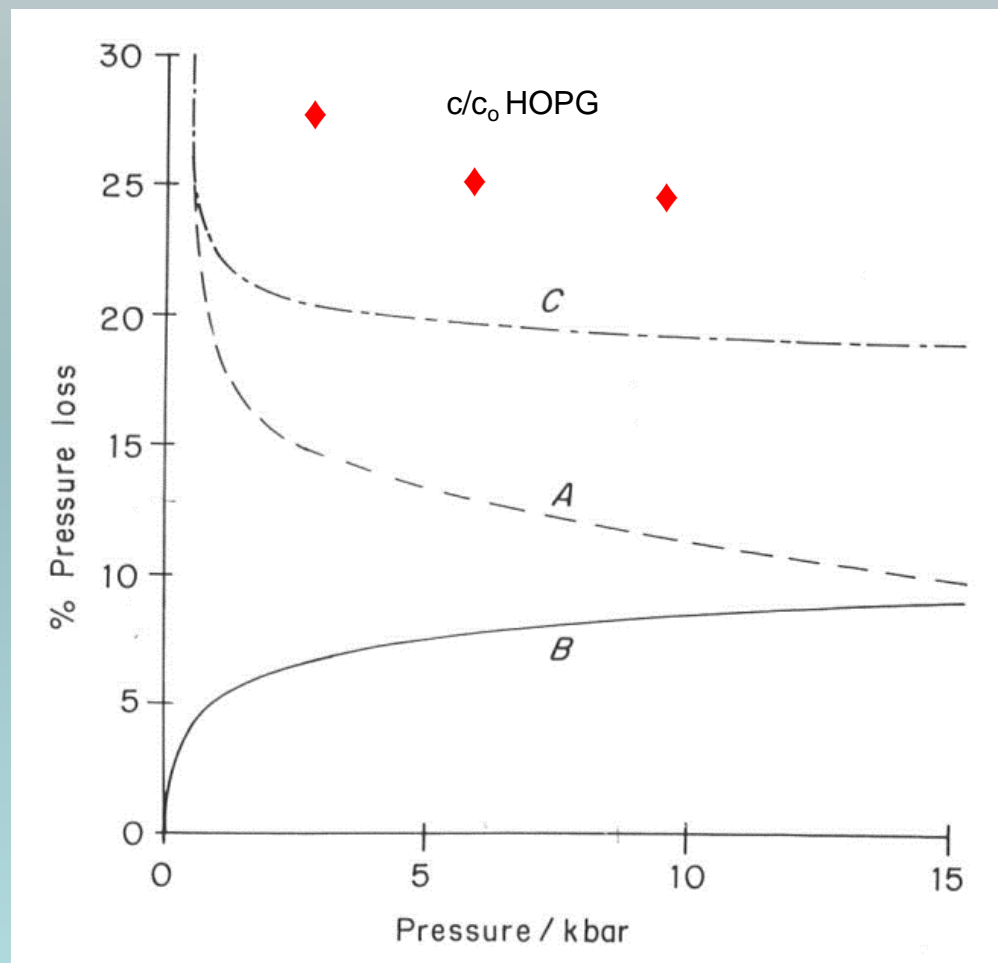

W.F. Sherman, A.A. Stadtmuller. *Exp. Tech. in H. Pressure Research.*

Results

Freezing and Cooling (VP_{constant}) $\rightarrow \Delta P \sim 25\%$

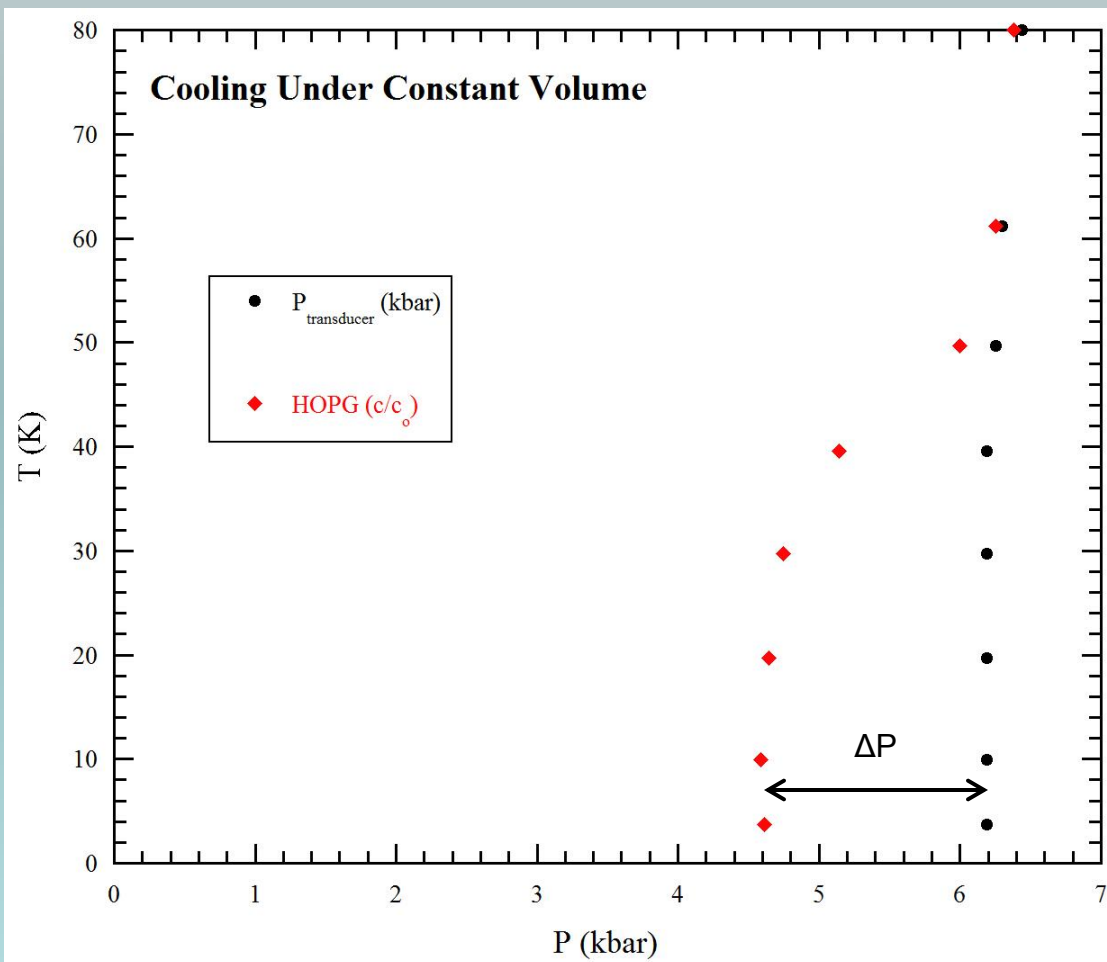


Results

Cooling to Freezing Point (P_{constant})

$\Delta P_{\text{AV}} \sim 5\%$

When systematically ensuring
that the pressure vessel is
completely full of solid He

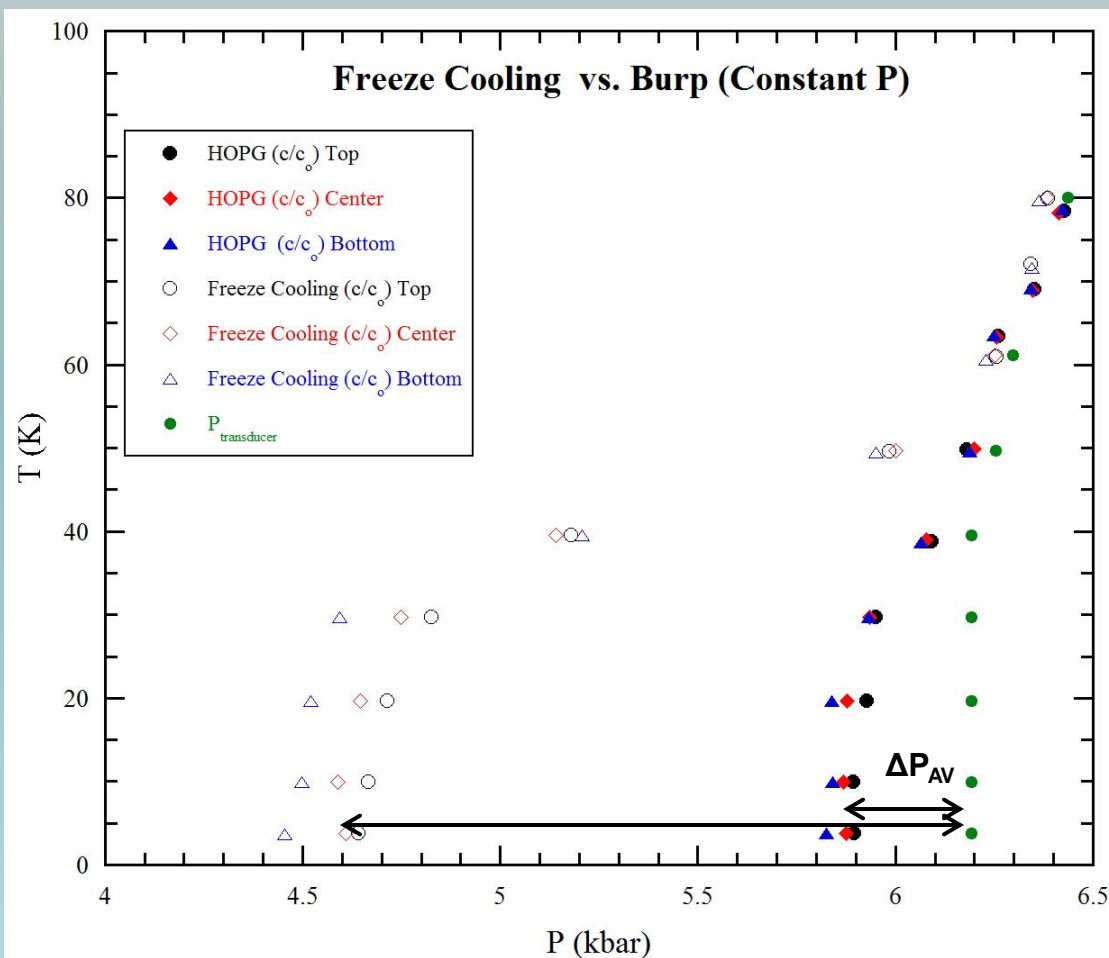


Results


Cooling to Freezing Point (P_{constant})

$$\Delta P_{AV} \sim 5\%$$

When systematically ensuring
that the pressure vessel is
completely full of solid He

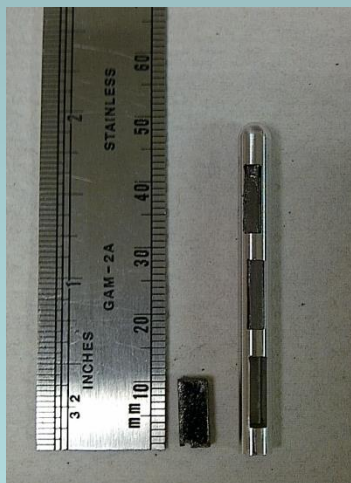

5mm x 10mm HOPG Xtals $\sim 10^0$ offset

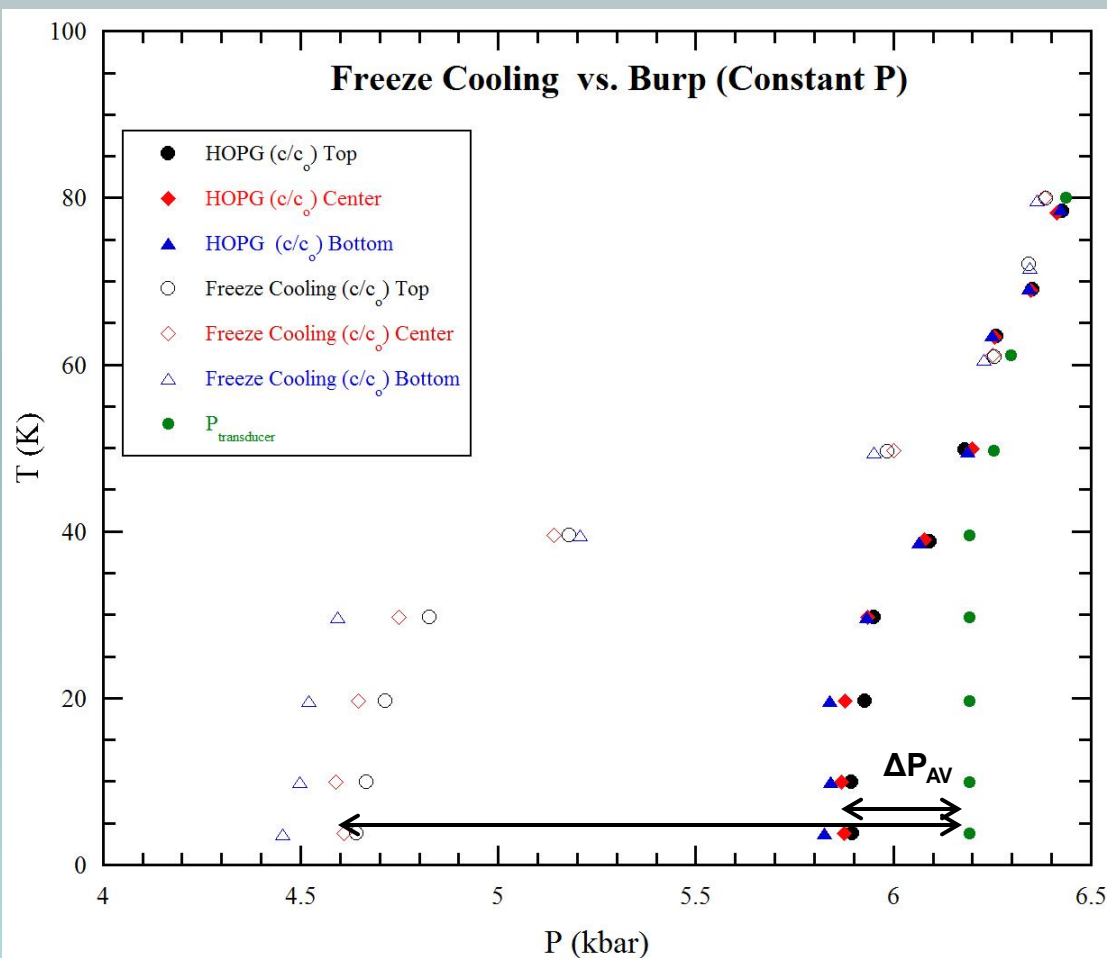

Results


◆ Freezing and Cooling (VP_{constant})

→ $\Delta P \sim 25\%$

● Freezing and Cooling (P_{constant}) for neutron elastic measurements of HOPG (c/c_0)

→ $\Delta P_{AV} \sim 5.5\%$

- A. Freezing (V_{constant})
- B. Cooling to Freezing Point (P_{constant})
- C. Freezing and Cooling (VP_{constant})

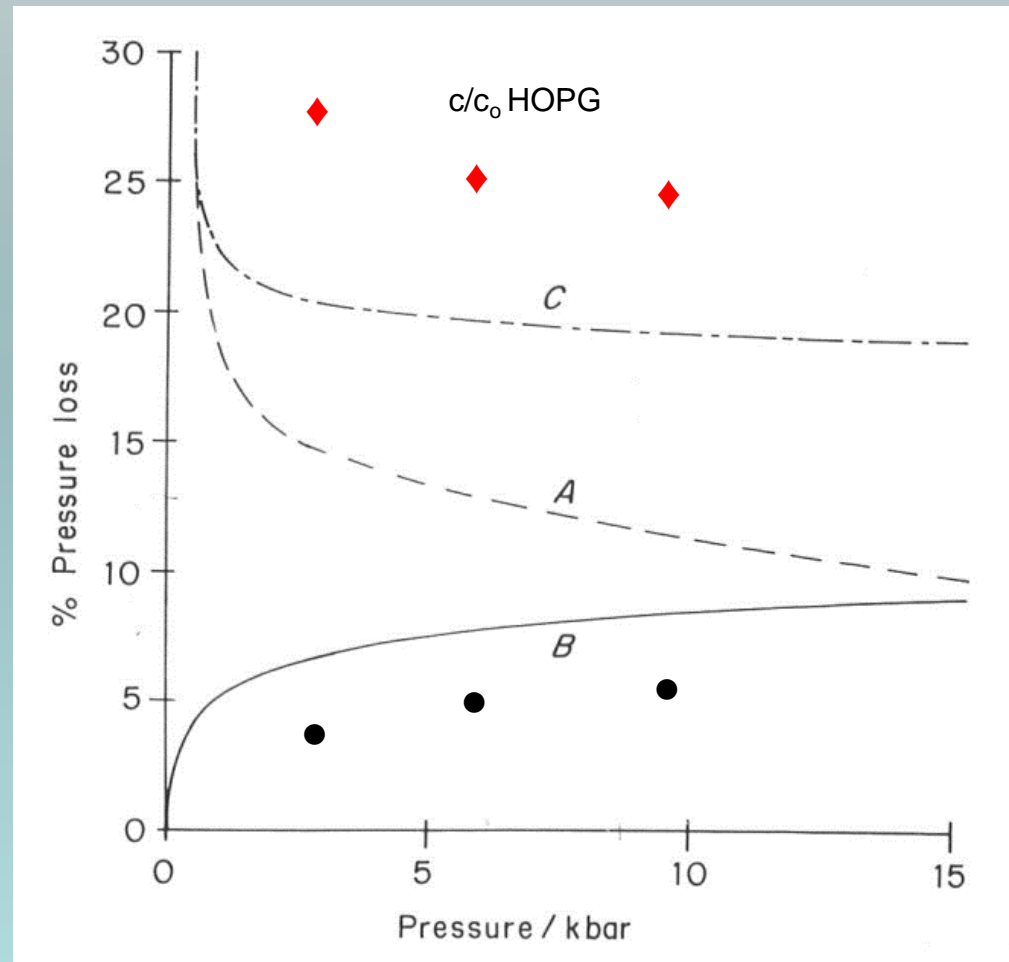


W.F. Sherman, A.A. Stadtmuller. *Exp. Tech. in H. Pressure Research.*

Results

Freezing and Cooling (P_{constant})
for neutron elastic measurements of
HOPG (c/c_0)

Freezing and Cooling (P_{constant}) $\rightarrow \Delta P_{\text{AV}} \sim 5\%$

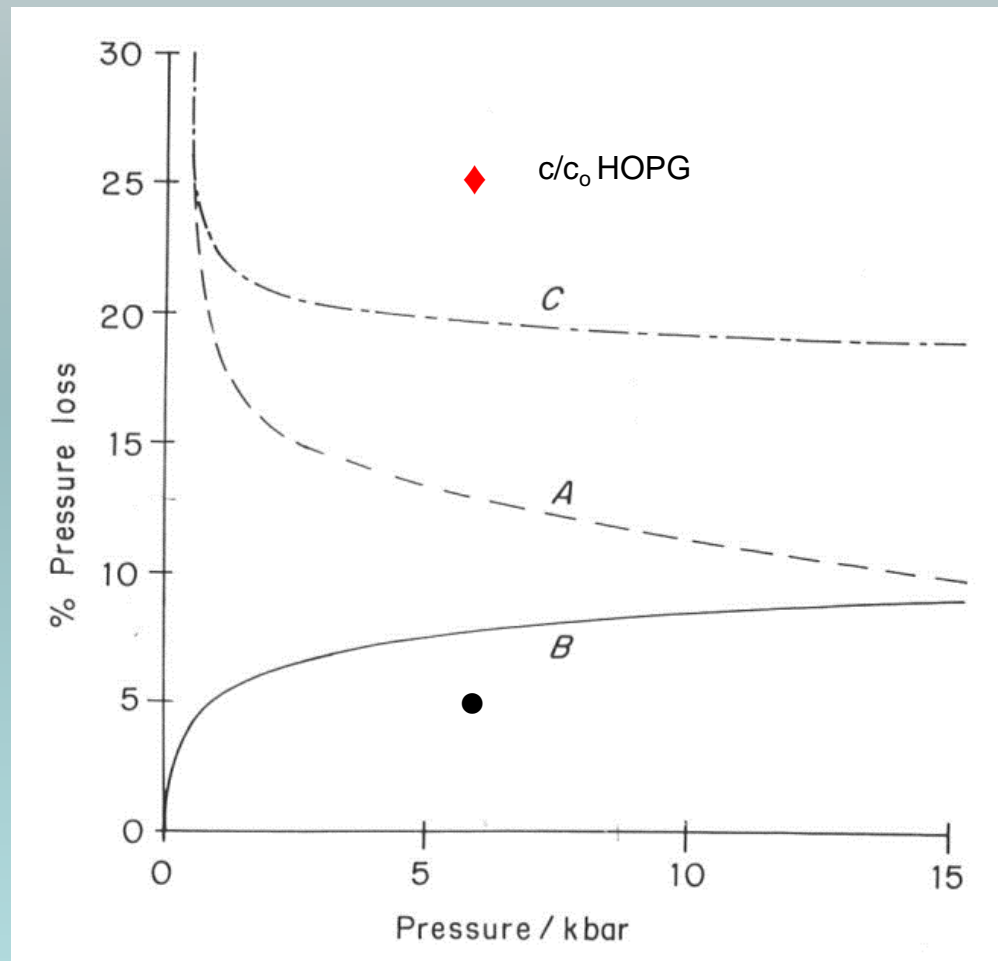


W.F. Sherman, A.A. Stadtmuller. *Exp. Tech. in H. Pressure Research.*

Results

Freezing and Cooling (P_{constant}) $\rightarrow \Delta P \sim 5\%$


When systematically ensuring
that the pressure vessel is
completely full of solid He

Bonus:
Sample pressure in-homogeneities
are minimized

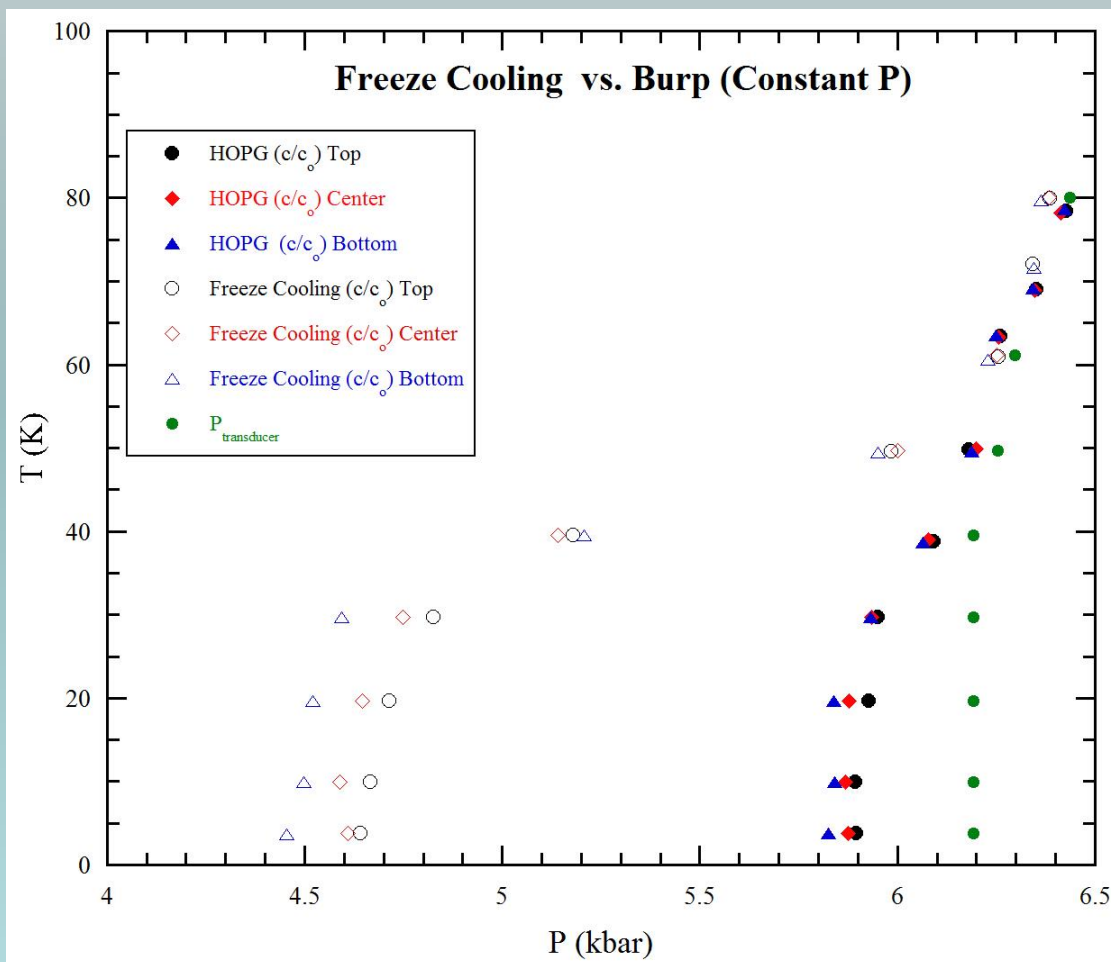


Results

Sample Pressure in-homogeneities (P_i) Comparison

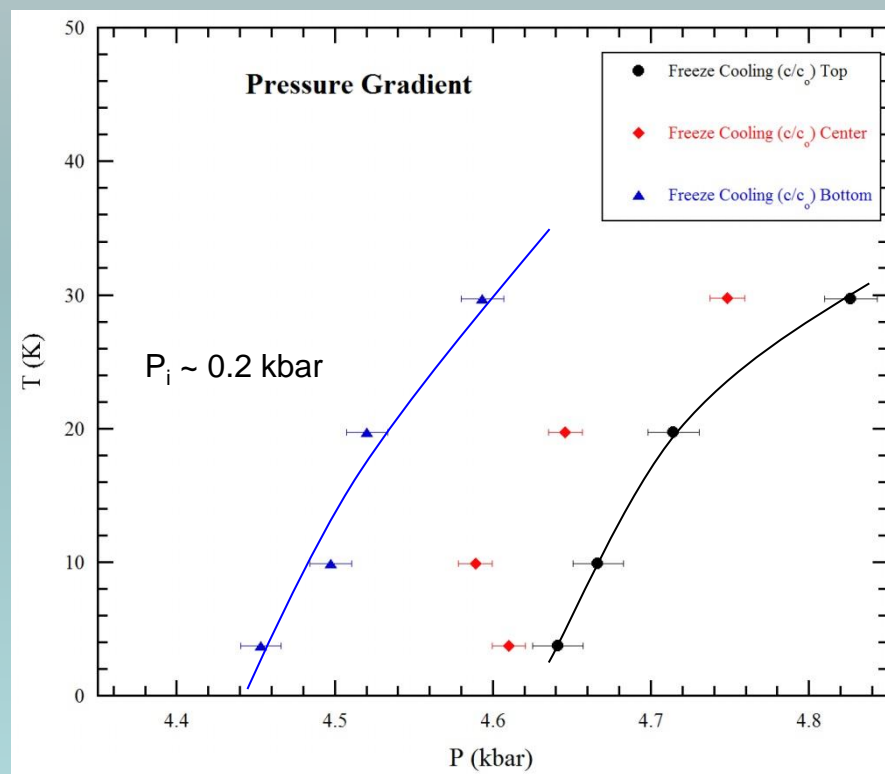

Freezing and Cooling (VP_{constant}) $\rightarrow \Delta P \sim 25\%$

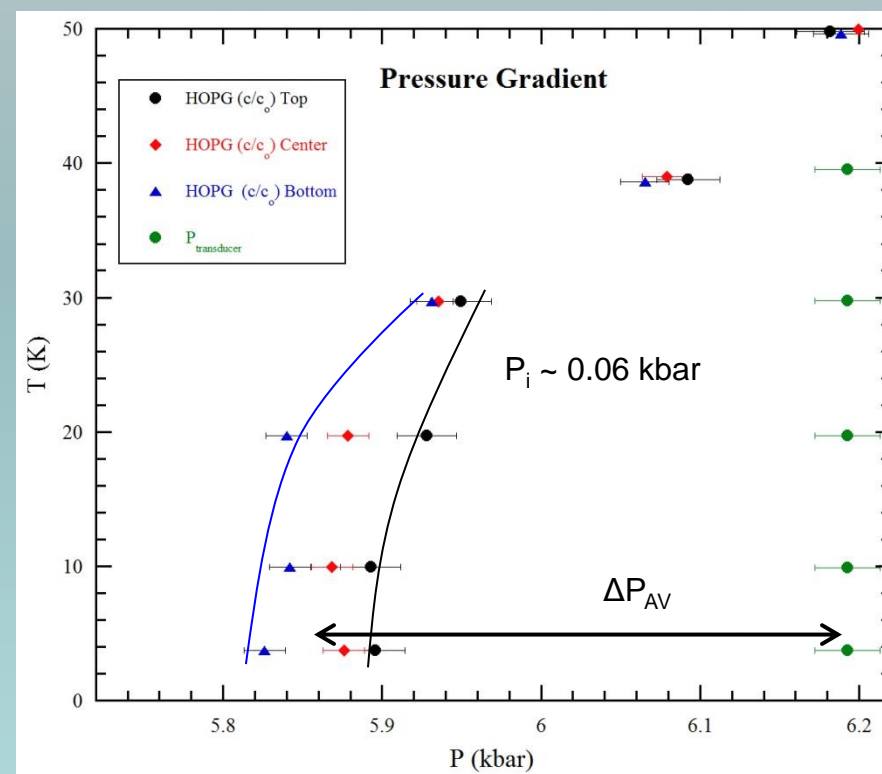



Burp Freezing and Cooling (P_{constant}) $\rightarrow \Delta P \sim 5\%$

Technique

- Use He as pressure media
- Monitor pressure line temperature
- Control line temperature well above desired PxT
- Cool down to a few Kelvin of PxT curve
- Apply pressure while heating line
- Begin ramp-cooling through the PxT curve
- Systematically heat line upon transducer pressure “freeze” providing enough power to counteract the cooling of the cell as noted in the sample stick sensor
- Once transducer pressure “melts” reduce line heater power to once again “freeze” the transducer reading
- Repeat until $P_{\text{Freeze}} \sim P_{\text{Melt}}$
- Continue cooling to base while still heating line
- Begin reducing line heater power when the cryostat cooling bottoms out
- Turn off line heater when $T_{\text{Line}} < T_{\text{Freeze}}$



Technique

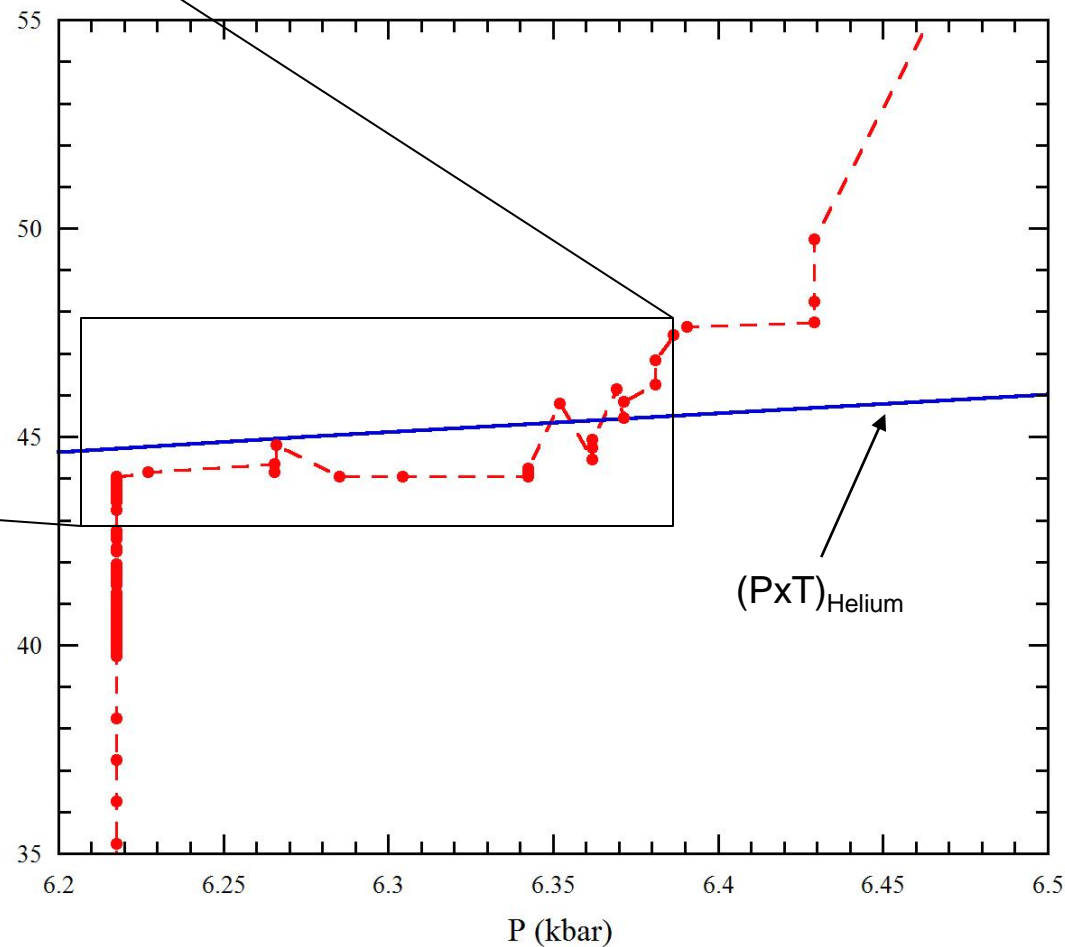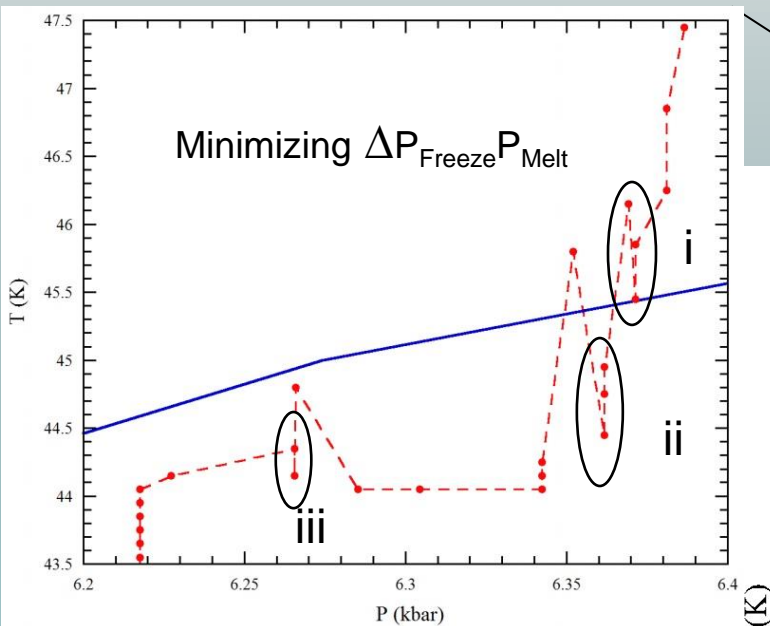


$$\Delta P_F P_M(i) < \Delta P_F P_M(ii) < \Delta P_F P_M(iii)$$

Results

Freezing and Cooling (P_{constant}) $\rightarrow \Delta P \sim 5\%$


When systematically ensuring that the pressure vessel is completely full of solid He

Bonus:
Sample pressure in-homogeneities are minimized

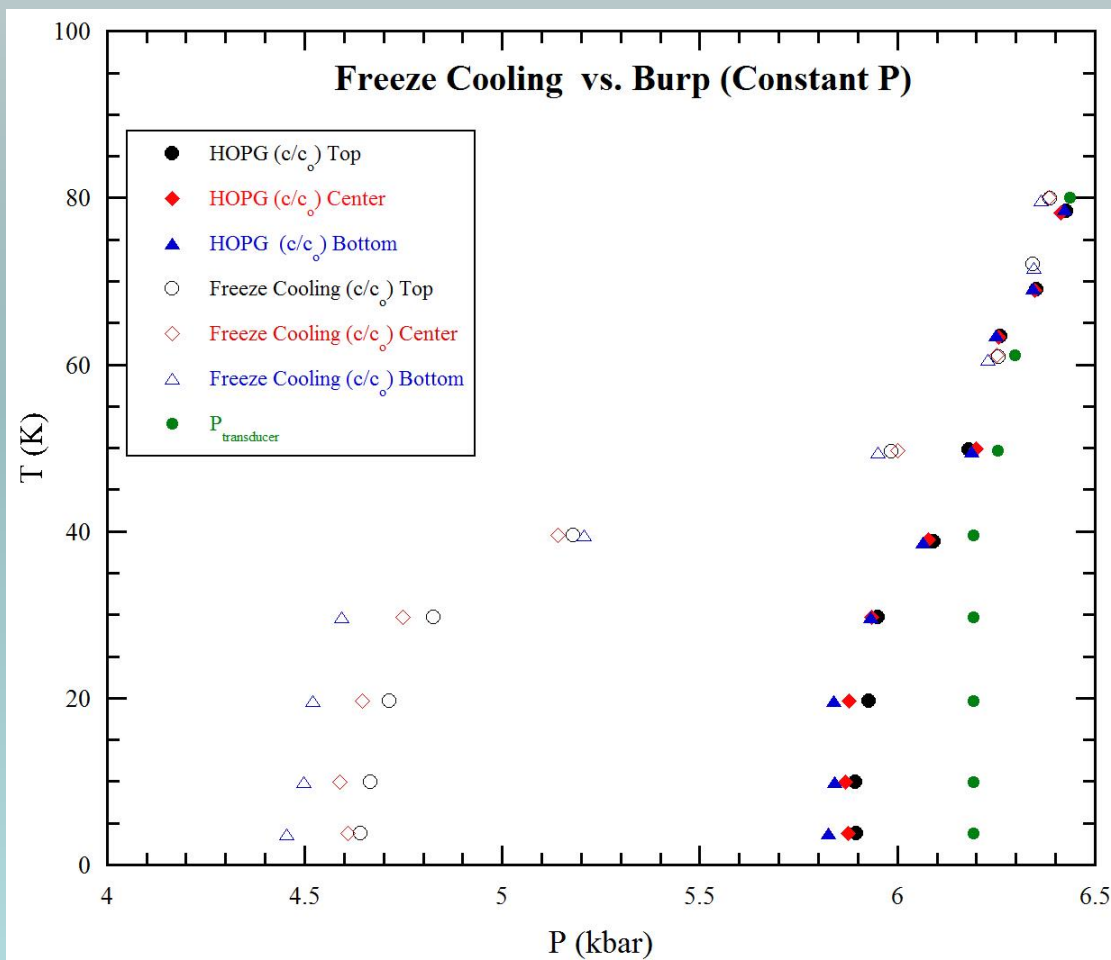


Results

Sample Pressure In-homogeneities (P_i) Comparison

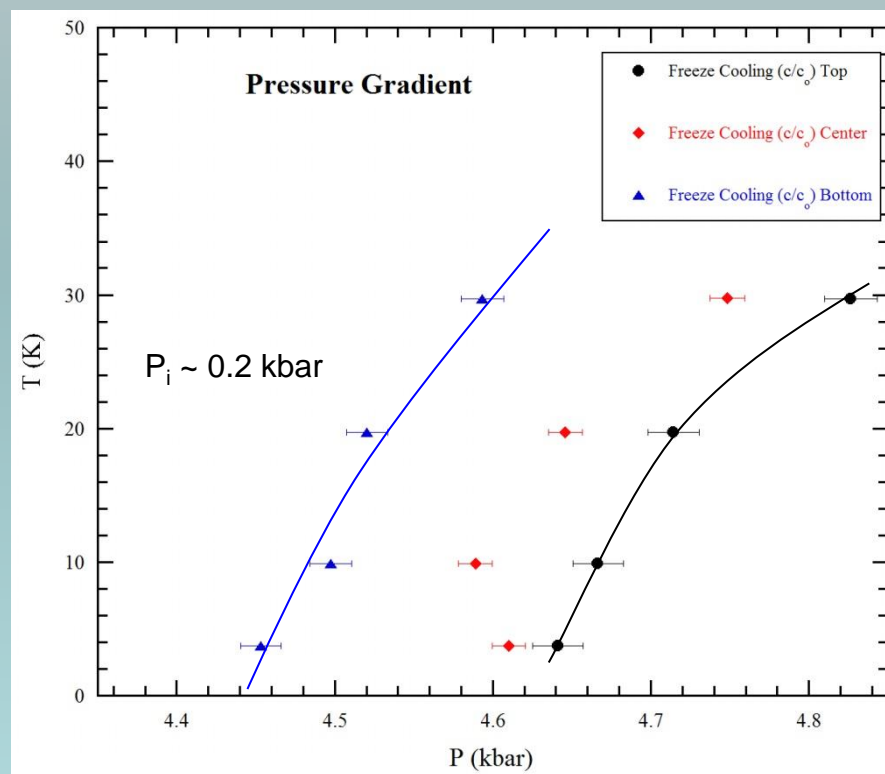

Freezing and Cooling (P_{constant}) $\rightarrow \Delta P \sim 25\%$

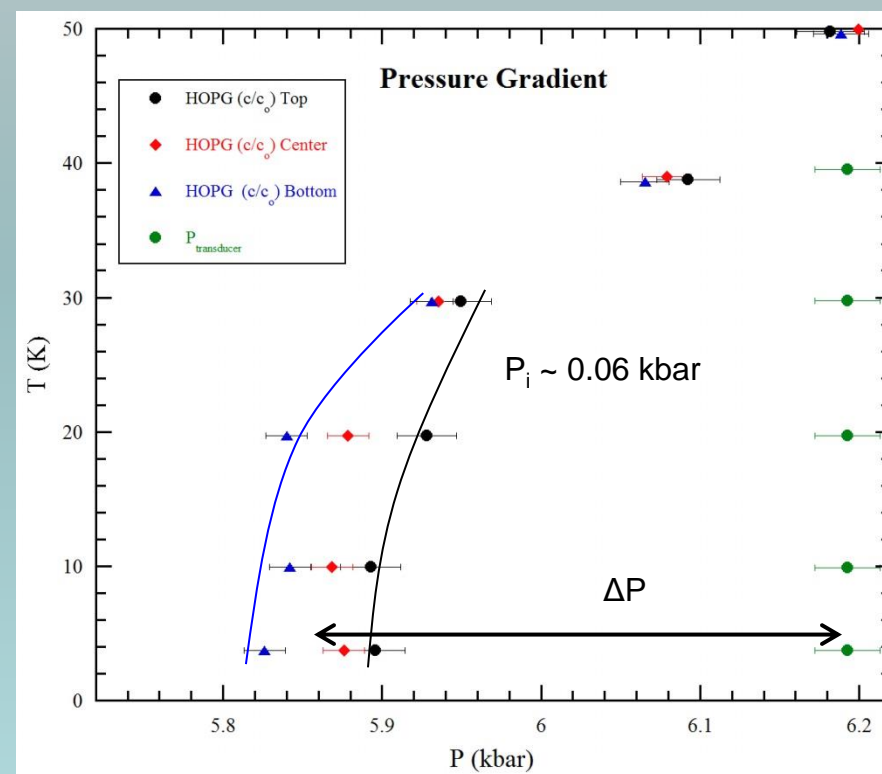

Burp Freezing and Cooling (P_{constant}) $\rightarrow \Delta P \sim 5\%$

Minimizing Pressure In-homogeneities for Large Samples in High Pressure Neutron Scattering Measurements

- ✓ Reasoning
- ✓ Apparatus
- ✓ Results
- ✓ Technique

Thank You!